


### LDS™ stapler:

The ligating and dividing system stapler is designed to ligate and divide blood vessels and is a very efficient and time saving device. The vessels are inserted in the end of the instrument jaws and on firing, two staples are inserted adjacent to each other while a knife simultaneously cuts the middle. The LDS™ stapler is of particular importance for ligation of mesenteric vessels and facilitating the removal of the spleen. This instrument is available as a re-usable, stainless steel model and as a totally disposable “Powered” LDS™. The disposable LDS™ is gas powered for ease of use. The disposable loading unit for the stainless steel LDS™ contains 15 pairs of staples.

### Main applications in veterinary surgery:

- Resection of omentum or omental adhesions
- Rapid ligation of vascular pedicles in a variety of procedures but especially useful for splenectomy

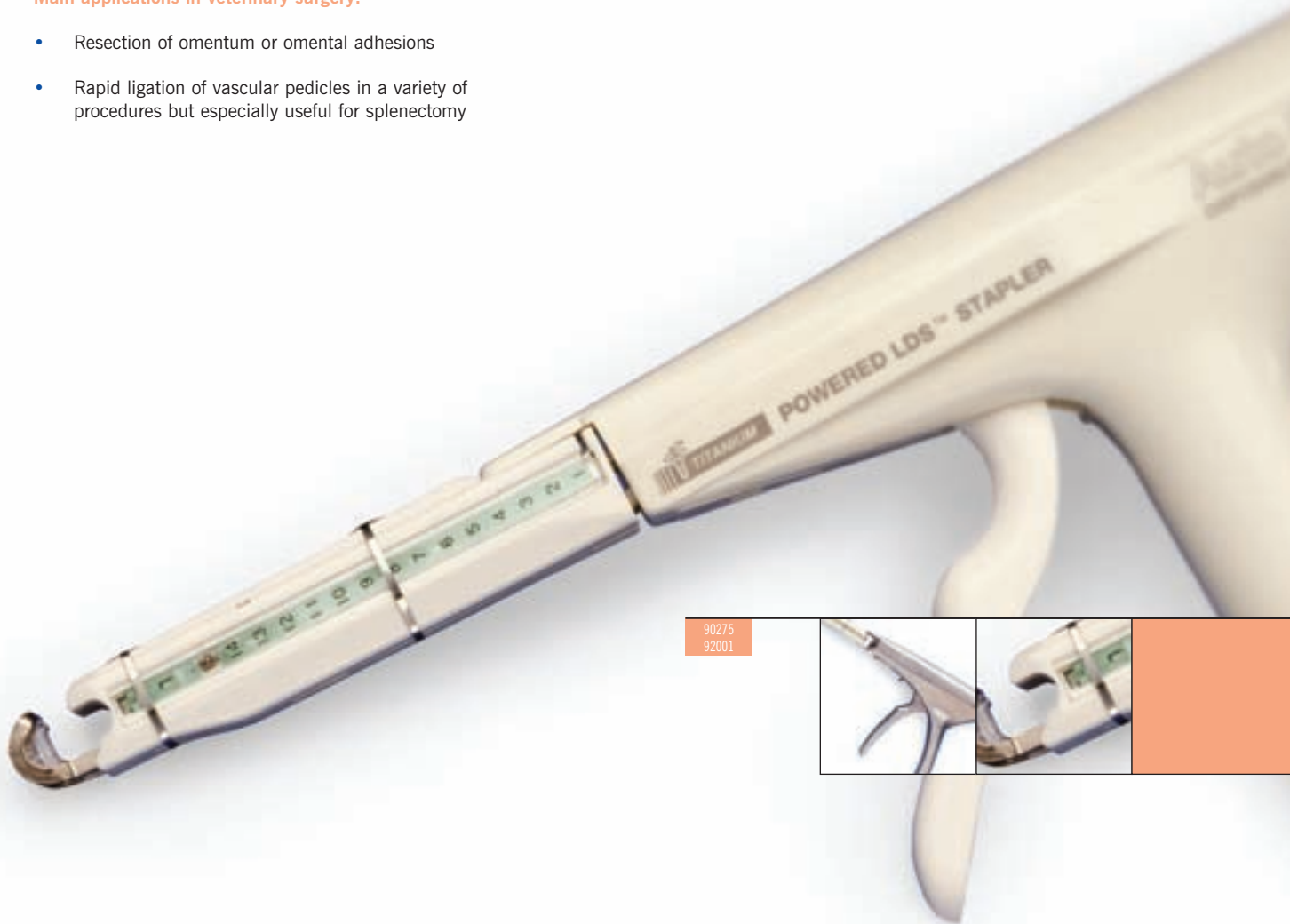

# LDS™ Staplers

# Premium Surgiclip™



The Premium Surgiclip™ contains 30 titanium superinterlocking clips which can be fired individually to ligate small blood vessels. Its design facilitates access to awkward areas. The superinterlock clip has a locking device which guarantees maximum security. This stapler has proved to be beneficial with exotic species where the vessels are small and access is often particularly difficult.